

LED BP MONITOR MF-09

MERCURY FREE BLOOD PRESSURE MONITOR

Measure
blood pressure,
to control
blood pressure.



Mercury is highly toxic heavy metal which if inhaled causes sleep disorder, agitation, paralysis and neurological disorder.

Monitor Blood Pressure with Medtech's Mercury Free LED BP Monitor.

KEY FEATURES

-  Battery 2 X AA
-  LED Display
-  USB port
-  Digital Sensor Technology

BENEFITS

-  Non Hazardous (Mercury free)
-  Accurate Reading
-  Light Weight

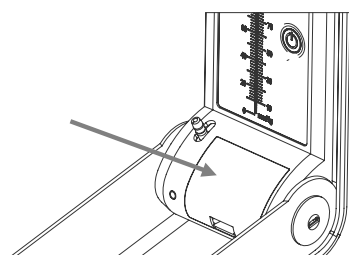
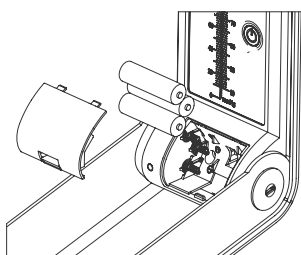
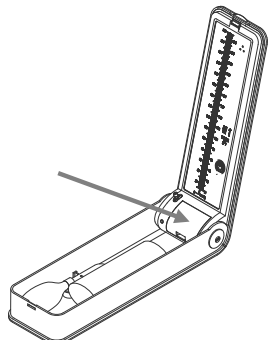


Useful For -

Measuring blood pressure to avoid cardiovascular disease, which is a leading cause of death.

HOW TO USE

Installation of Batteries:

| | | |
|--|--|---|
| <p>1. Remove the battery cover.</p>  | <p>2. Insert a new set of batteries into the battery compartment according to the polarities (+) and (-) are correct. Use only LR6, AA or equivalent batteries bag.</p>  | <p>3. Close the battery cover.</p>  |
|--|--|---|

4. Changing Battery pack: When battery low indication starts blinking continuously in Red color, change it with new battery pack in time.
 5. Remove batteries when device is not in use for longer period of time.

Specification

| | |
|----------------------------------|---|
| Measurement method | Auscultation using a stethoscope |
| Pressurization | Manually by Rubber bulb |
| Accuracy | Pressure +/- 3mmHg |
| LED column display Resolution | 2mmHg |
| Power supply | 3 x Alkaline battery (type LR6, type AA) or Equivalent or recommended DC supply |
| Deflation | Using control valve |
| Operating Temperature conditions | +10°C to +40°C |

Please take advice of the authorized medical practitioner before use.

Note: The information provided is best to our knowledge and subject to change without notice. Picture for illustration purposes only, actual product and accessories may vary. Specifications are subject to change without notice.

OUR PRODUCT RANGE

Nebulizer | BP Monitor (Digital/LED) | Steam Inhaler | | Thermometer (Standard / Infrared) | Pulse Oxymeter
 Oxygen Concentrator | Air Pressure Mattress | Auto CPAP | Needle Burner & Syringe Destroyer

MEDTECH LIFE PVT. LTD. (Formerly Nulife Global Medical Devices Pvt. Ltd.)

B6, Byculla Service Industries, D.K. Marg, Byculla, Mumbai 400 027, India.

+9122 2374 8372 / 73 +9122 6657 8989 info@medtechlife.com

Domestic: +91 82917 26885/ +91 82917 26884 / 83

Export: +91 91369 23382 / +91 87676 84540 / +91 93263 56572